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8th inter-laboratory study on the detection of Verocytotoxin-producing *E. coli* (VTEC) belonging to the serogroups most involved in human infections (O157, O111, O26, O103, and O145) in water samples

Laboratory procedure

Introduction

In the absence of specific international standards for the detection of VTEC in water, the method developed for the detection of VTEC in foodstuffs, ISO/TS 13136 “Microbiology of food and animal feed — Real-time polymerase chain reaction (PCR)-based method for the detection of food-borne pathogens — Horizontal method for the detection of Shiga toxin-producing *Escherichia coli* (STEC) belonging to O157, O111, O26, O103 and O145 serogroups”, has been adapted to this particular matrix.

The procedure comprises the following sequential steps:

- Filtration of water samples through membranes of mean pore size 0.45 µm
- Transfer of the membranes to the enrichment medium
- Microbial enrichment
- Nucleic acid extraction
- Detection of virulence genes
- Detection of serogroup-associated genes
- Isolation from positive samples.

1. Filtration of the water samples and transfer of the membranes to the enrichment medium

These steps are performed according to ISO 8199, section 8.2.5. Briefly, a sterile filtration apparatus connected to a source of vacuum is used. Place a sterile membrane filter of mean pore size 0.45 µm on the porous disc of the filter base, with only the outer part of the membrane filter being grasped by flat-ended sterile forceps. Pour the content of the bottle containing the test sample and apply sufficient vacuum (about 70 kPa) to filter the water

through the membrane. The membrane is then placed in a flask containing 100 ml of the enrichment liquid culture medium. This is constituted by buffered peptone water (BPW), which is indicated to analyze samples that are supposed to contain stressed target bacteria. The flask is then incubated for 18 -24 h at 37 °C ± 1°C

2. Nucleic acid extraction, detection of virulence and serogroup-associated genes, and isolation of VTEC

These steps are performed according to the method CEN/ISO TS 13136 “Microbiology of food and animal feed — Real-time polymerase chain reaction (PCR)-based method for the detection of food-borne pathogens — Horizontal method for the detection of Shiga toxin-producing *Escherichia coli* (STEC) belonging to O157, O111, O26, O103 and O145 serogroups”.

This method is described in the **Annex**. Briefly, it is based on a Real-time PCR screening of enrichment cultures to detect the presence of virulence genes (*vtx1* and *vtx2*, and *eae*) and serogroup-specific genes for O26, O103, O111 and O145, followed by the isolation of VTEC from PCR-positive samples, accomplished by an immuno-concentration enrichment step specific for the serogroups identified in the PCR step. One ml aliquot of the enrichment culture performed as described above is used for DNA extraction and purification, accomplished according to the ISO 20837 “Microbiology of food and animal feeding stuffs - Polymerase chain reaction (PCR) for the detection of foodborne pathogens - Requirements for sample preparation for qualitative detection”. The remaining culture shall be stored at 4°C for the isolation steps that will follow a positive PCR result. The method is sequential:

Step 1: Detection of the genes *vtx1*, *vtx2* and *eae*.

Step 2: Samples positive for both *vtx* and *eae* are tested for the serogroup-associated genes (molecular serogrouping).

Step 3: Isolation of the VTEC strain; samples positive at the same time for *vtx*, *eae*, and at least one of the serogroup-associated genes are submitted to a further step aimed at isolation of the VTEC strain. This requires serogroup-specific enrichment based on IMS or other immuno-capture suitable approaches. A guideline for the isolation of the different VTEC serogroups is also included in **Annex 1**.

Step 4: Characterization of the isolate i.e. identification, detection of *vtx* genes, the *eae* gene and the serogroup gene.

Annex 1

EU Reference Laboratory for *E. coli*
Department of Veterinary Public Health and Food Safety
Unit of Foodborne Zoonoses
Istituto Superiore di Sanità

Laboratory procedure for the detection of Shiga toxin (Verocytotoxin)-producing *Escherichia coli* (STEC/VTEC) belonging to O157, O111, O26, O103 and O145 serogroups - Qualitative Method

Introduction

Shiga toxin-producing *Escherichia coli* (STEC) cause severe disease in humans such as haemorrhagic colitis and haemolytic uraemic syndrome (HUS). Although STEC may belong to a large number of serogroups, those that have been firmly associated with severe human disease, in particular HUS, belong to O157, O26, O111, O103, O145 (1), and represent the targets of this Laboratory Procedure.

In this Laboratory Procedure, the wording Shiga toxin (Stx) is synonymous of Verocytotoxin (Vtx).

The following nomenclature has been adopted throughout the text:

stx: Shiga toxin genes (synonymous of *vtx*)

Stx: Shiga toxin (synonymous of *vtx*)

STEC: Shiga toxin-producing *Escherichia coli* (synonymous of VTEC: Verocytotoxin-producing *Escherichia coli*).

1. Scope

This Laboratory Procedure describes a horizontal method for the detection of (i) the major virulence genes of STEC (2,3), and (ii) the genes associated with the serogroups O157, O111, O26, O103 and O145 (3,4).

In the case of detection of these genes, the isolation of the strain is attempted, to confirm the simultaneous presence of the genes in the same live bacterial cell.

The analytical approach is based on the use of the Real-time PCR..

This Laboratory Procedure is applicable to any enrichment culture from:

- products intended for human consumption and the feeding of animals,
- environmental samples in the area of food production and food handling,
- environmental samples in the area of primary production.
- Water samples

2. Normative references

ISO/DIS 7218

General requirements and guidance for microbiological examinations

ISO/DIS 20837

Requirements for sample preparation for qualitative detection

ISO/DIS 20838

Requirements for amplification and detection of qualitative methods

ISO/DIS 22174

General method specific requirements

3. Terms and definitions

3.1 Shiga toxin-producing *Escherichia coli* (STEC): Microorganism possessing the Stx-coding genes

3.2 Shiga toxin-producing *Escherichia coli* (STEC) potentially pathogenic to humans: Microorganism possessing the Stx-coding genes and the intimin-coding gene *eae*.

3.3 Shiga toxin-producing *Escherichia coli* (STEC) highly pathogenic to humans: Microorganism possessing the Stx-coding genes, the intimin-coding gene *eae* and belonging to one of the serogroups in the scope of the present Laboratory Procedure.

4. Principle

4.1 General

The detection of STEC and of the 5 serogroups comprises the following sequential steps:

1. Microbial enrichment

2. Nucleic acid extraction
3. Detection of virulence genes
4. Detection of serogroup-associated genes in samples positive to point 3
5. Isolation from samples positive to points 3 and 4.

A flow-diagram of the whole procedure is given in Annex A.

4.1.1 Microbial enrichment

The number of STEC cells to be detected is increased by incubating the test portion in a non-selective liquid nutrient medium, either:

- (a) Tryptone-soy broth (TSB) supplemented with 1.5 gr/l bile salts n. 3 and 16 mg/l of novobiocin (mTSB+N).
- (b) Buffered peptone water (BPW)
- (c) Tryptone-soy broth (TSB) supplemented with 1.5 gr/l bile salts n. 3 and 12 mg/l of acriflavin (mTSB+A) for analysis of dairy products.

The broth mTSB is to be used when analysing matrices suspected to contain high levels of contaminating microflora. Novobiocin and acriflavin inhibit the growth of Gram-positive bacteria and promote the growth of Gram-negative cells including STEC. The broth BPW is to be used to analyse samples which are supposed to contain stressed target bacteria (such as frozen products), to resuscitate stressed STEC cells, and expected lower levels of contaminating microflora than in fresh samples.

NOTE: The addition of novobiocin is controversial and has been investigated by several authors. It has been observed that the Minimum Inhibitory Concentration (MIC) of the antibiotic for non-O157 STEC is lower than for O157 strains (5). The addition of novobiocin in the enrichment broth mTSB at the usual concentration of 20 mg/l, as specified in the ISO 16654:2001 standard, seems to inhibit the growth of about one third of non-O157 strains (6) increasing the risk of false negative results. The 16 mg/l concentration represents the most balanced option to inhibit the growth of the contaminating background microflora expected in enrichment cultures of food samples, while still allowing the growth of STEC cells.

4.1.2 Nucleic acid extraction

Bacterial cells are separated from the enrichment medium and lysed. The nucleic acid is then purified according to any protocol allowing the production of template DNA suitable for Real Time PCR.

4.1.3 Target genes

The purified nucleic acid is used for the detection of the following target genes:

- the main virulence genes for STEC: *stx* genes, encoding the Shiga toxins and the *eae* gene, encoding a 90KDa protein, the intimin, involved in the attaching and effacing mechanism of adhesion, a typical feature of the pathogenic STEC strains.
- the *rfbE* (O157), *wbdI*(O111), *wzx*(O26), *ihp1*(O145) and *wzx*(O103) genes, to identify the corresponding serogroups.

4.1.4 Detection

The Real-time PCR products are detected by light emission in a 5' nuclease PCR assay.

4.1.5 Isolation

Once a STEC is detected and the serogroup identified, in order to isolate the STEC strain, a serogroup-specific enrichment is performed followed by plating onto the agar Tryptone-bile-glucuronic medium (TBX) or onto a specific selective medium where available (see note in Annex A).

5. Diluent, culture media and reagents

5.1 Culture media

5.1.1 modified Tryptone- Soy broth (mTSB)

5.1.1.1 Basic medium

Composition

Casein peptone	17 g
Soy peptone	3 g
D(+) Glucose	2.5 g
Sodium chloride	5 g
Dipotassium hydrogen phosphate (K ₂ HPO ₄)	4 g
Bile salts no. 3	1.5 g

Distilled water to 1 L

pH 7.4 ± 0.2

Preparation

Dissolve the components or the dehydrated medium in water. Adjust pH with a pH-meter pH 7.4 +/- 0.2 at 25°C (6.8 and sterilize by autoclaving at 121°C 15 min (6.9).

5.1.1.1.2 Novobiocin solution

Composition

Novobiocin	0,16	g
Water	10	ml

Preparation

Dissolve the novobiocin in the water and sterilize by membrane filtration.

Prepare on the day of use.

5.1.1.1.3 Acriflavin solution

Composition

Acriflavin	0,12	g
Water	10	ml

Preparation

Dissolve the acriflavin in the water and sterilize by membrane filtration.

Prepare on the day of use.

5.1.1.1.4 Preparation of the complete medium

Immediately before use, add 1 ml of novobiocin (5.1.1.1.2) or acriflavin solution (5.1.1.1.3) to 1000 ml of cooled mTSB (5.1.1.1.1)

The final concentration of novobiocin is 16 mg per litre of mTSB.

The final concentration of Acriflavin is 12 mg per litre of mTSB

5.1.2 Buffered peptone water (BPW)

Peptone	10	g
Sodium chloride	5.0	g
Disodium phosphate (Na ₂ HPO ₄)	3.5	g
Potassium dihydrogen phosphate (KH ₂ PO ₄)	1.5	g
Water	to 1000	ml

pH 7.2 ± 0.2

Preparation

Dissolve the components or the dehydrated powder in the water. Adjust pH with a pH-meter pH 7.4 ± 0.2 at 25 °C (6.8) and sterilize by autoclaving at 121 °C 15 min (6.9).

5.2 Reagents for nucleic acid extraction

The reagents to be used for nucleic acid extraction are not listed being dependent on the method adopted (9.3)

5.3 Reagents for PCR

See ISO 20838.

5.3.1 Oligonucleotides (primers) and detection probes

Primers and probes for specific detection of the target gene sequences by Real time PCR are listed in Annex E.

6. Equipment

Usual microbiological laboratory equipment (see ISO 7218) and, in particular, the following:

6.1 Water bath up to 100 °C

6.2 Incubator according to ISO 7218, 37 ± 1 °C

6.3 Apparatus for nucleic acid extraction

Appropriate equipment according to the method adopted. If needed.

6.4 Pipettes for volumes between 1 µl and 1000 µl

6.5 Thin walled Real-Time PCR microtubes (0,2 ml /0,5 ml reaction tubes), multi-well PCR microplates or other suitable light transparent disposable plasticware.

6.6 Thermal cycler Several brands of apparatus are available and can be chosen according to the laboratory policies.

6.7 Apparatus for detection of the PCR product

Light emission following 5' nuclease PCR assay is detected by the Real-time PCR apparatus

6.8 pH-meter capable of measuring to an accuracy of ± 0.05 pH units and its resolution shall be 0.01 pH units

6.9 Autoclave according to ISO 7218

6.10 Stomacher peristaltic blender with sterile bags possibly with a device for adjusting speed and time

7. Sampling

It is important that the laboratory receives a sample that is truly representative and has not been damaged or changed during transport or storage. Sampling is not part of the method specified in this Technical Specification. See the specific International Standard dealing with the product concerned. If there is no specific International Standard dealing with sampling of the product concerned, it is recommended that the parties concerned come to an agreement on this subject.

8. Preparation of test sample

Prepare the test sample in accordance with the specific International Standard dealing with the product concerned where available, or other suitable guidelines.

9. Procedure

9.1 Test portion and initial suspension

9.1.1 General

Use the necessary quantity of enrichment medium to yield a final dilution of 1/10 of the original test portion. Filters used for membrane filtration of water samples are added to 100 ml of enrichment broth (BPW).

9.1.2 For matrix sample supposed to contain high level of annex floras

For solid matrices, aseptically transfer a test portion of sample to a stomacher bag containing the appropriate amount of mTSB added with novobiocin or acriflavin (5.1.1.1.4). Bags with filters should be preferred.

Homogenise in a stomacher (see ISO 7218) (6.10).

For liquid matrices, transfer the test portion of liquid sample, using a sterile pipette, directly into the base of the tube/bottle containing the enrichment broth mTSB added with novobiocin or acriflavin (5.1.1.1.4).

9.1.3 For matrix sample supposed to contain stressed target bacteria

In the case of frozen products, allow them to thaw at room temperature, then transfer the test portion to a stomacher bag containing the appropriate volume of BPW (5.1.2) and proceed as above.

9.2 Enrichment

9.2.1 Incubation

Incubate the stomacher bag or the tube/bottle (9.1.2) at $37\text{ °C} \pm 1\text{ °C}$ for 16 h to 18 h.

9.2.2 Process control (for Real-time PCR)

Perform process control according to ISO 22174.

Guidance on internal amplification control (IAC) and process control are given in Annex D.

9.3 Nucleic acid preparation

Use an appropriate nucleic acid extraction procedure for Gram-negative bacteria. A comprehensive collection of methods can be found in (7). Alternatively, commercial kits may be used according to the manufacturers' instructions.

9.4 Real-time PCR amplification

9.4.1 General

Follow all requirements for the PCR amplification as described in ISO 20838 "Microbiology of food and animal feeding stuffs – Polymerase chain reaction for the detection of food pathogens – amplification and detection".

Primers and detection probes for conducting the real-time PCR are described in Annex E.

9.4.2 PCR controls

In accordance with ISO 22174, examples of PCR controls are given in Annex D.

9.4.3 Detection of PCR products

Light emission is captured by the apparatus once generated during the amplification.

9.4.4 Interpretation of PCR results

The PCR results obtained, including the controls specified in ISO 22174 and in Annex D, are interpreted by the software linked to the apparatus. During amplification, the Software monitors 5'-nuclease PCR amplification by analysing fluorescence emissions (Rn) of the reporter dye for each sample. ΔRn was Rn minus the baseline reporter dye intensity established in the first few cycles. At the end of the PCR, a reaction was considered positive if its ΔRn curve exceeded the threshold, defined as 10 times the standard deviation of the mean baseline emission calculated between the first few cycles. The cycle threshold (Ct) was defined as the cycle number at which a sample's ΔRn fluorescence crossed the determined threshold value.

If the results are ambiguous, check the emission curves. Positive samples give a curve with a clear increase in fluorescence, starting from a number of cycles corresponding to the Ct .

If the controls yield unexpected results, repeat the procedure.

9.5 Strain isolation

The pattern of virulence genes, exhibited by STEC which are considered to be pathogenic to humans, is complex and it is possible that different strains presenting only part of the virulence gene pattern may be present in the food sample at the same time. Therefore, the isolation of the STEC strains is required to confirm that the positive PCR signals are generated from genes that are simultaneously present in the same live bacterial cell.

A serogroup-specific enrichment followed by direct plating onto suitable solid media and screening of the colonies for the presence of the virulence genes are required.


The real-time PCR protocol described in this Laboratory Procedure may be used in order to confirm the presence of the virulence genes in the isolated colonies. In alternative other equivalent PCR protocols can be downloaded from the EU RL VTEC website (<http://www.iss.it/vtec/work/cont.php?id=152&lang=2&tipo=3>).

STEC isolation is described in the flow chart of **Annex B**.

NOTE: In the absence of positive results for the presence of *eae* gene, and generally in absence of positivity to serogroup-associated genes isolation is not attempted.

Annex A

Flow-diagram of the screening procedure

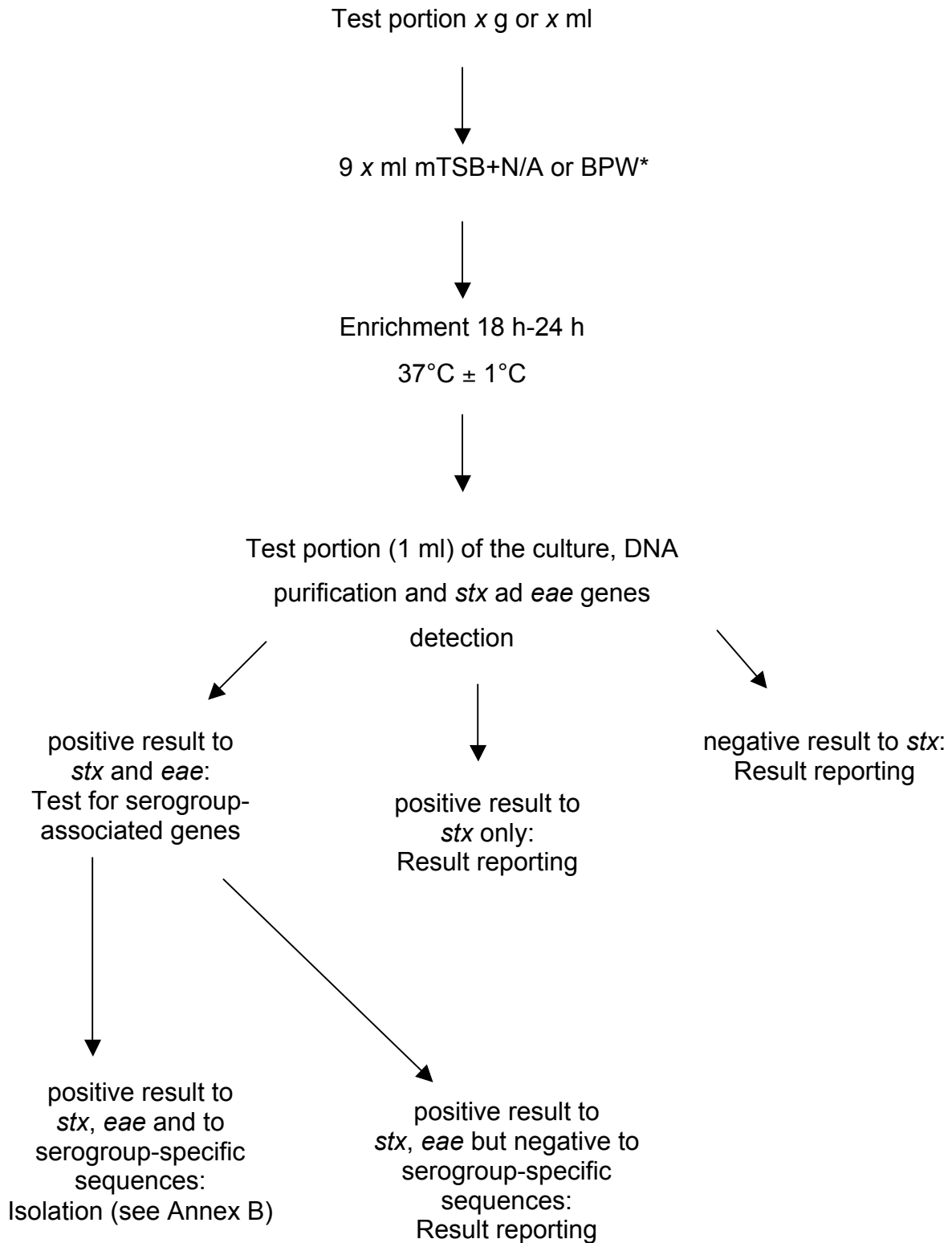
*= Not applicable to water samples processed by membrane filtration method

Annex B


Flow-diagram of the isolation procedure

Serogroup-Specific Enrichment (SSE)

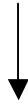Enrichment broth streaked onto suitable solid media.


Incubation for 18 h to 24 h at 37°C ± 1°C

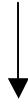
Pick up to 50 colonies with *E. coli* morphology. Point-inoculate on Nutrient Agar (NA) (single colonies) and H₂O (5 pools by 10 colonies each). Perform *stx* and *eae* detection on isolated colonies or pools.


If a colony is positive for the presence of the genes identified at the screening step go to the following step. If a pool is positive, incubate NA. Test individual colonies composing the positive pool as above.

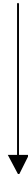Identify positive colonies as *E. coli* and check the serogroup (e.g. by PCR or slide agglutination)

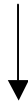Further characterisation (optional): send the strain to a Reference Laboratory

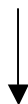Result reporting

Annex D Internal Amplification Control

Three different internal amplification controls (IACs) can alternatively be used in Real-time PCR:

- A commercially available TaqMan[®] Exogenous Internal Positive Control (Applied Biosystems, Foster City, CA, USA). The reagent kit include all reagent necessary (primers, a Vic[™]probe, IAC target DNA and blocking solution). The IAC target DNA must be diluted 10 times to achieve a copy number of approximately 100 per PCR reaction. The PCR product length is not declared to the customer.
- The open formula pUC 19 based internal amplification control IAC developed by Fricker et al. (11). Approximately 100 copies of target DNA (pUC 19) should be used per PCR reaction. The size of the IAC was 119 bp.
- A recombinant plasmid (named pIAC-STEC) can be used in the *stx*-specific real-time PCR assay (Auvray *et al.* personal communication). This IAC contains the following DNA fragment cloned into the *EcoRI* site of pUC19:5'-
ATTTTTGTTACTGTGACAGCTGAAGCTTTACGTGAATCGCCAGCGGCATCAGC
ACCTTGTCGCCTTGCGTATAGATGTTGATCTTACATTGAACTGGGGAATT-3'
(bold letters: *stx1/stx2* forward and reverse primers binding sites sequences; underlined sequence: IAC-probe binding site). The IAC is co-amplified with *stx* genes using the same primers as *stx* (Annex E), under the same conditions and in the same PCR tube (12). It is detected by a specific DNA probe (5' [Red640]-CAAGGCGACAAGGTGCTGATGCCG-[BHQ2] 3') included in the PCR mix at a concentration identical to that of the *stx1* and *stx2* DNA probes (both labelled with [Fluorescein] and [BHQ1] at their 5' and 3' ends, respectively). 64 copies of the IAC should be used per PCR reaction. The IAC PCR product is 96 bp-long. The performance of the resulting *stx*-IAC real-time PCR assay has been shown using artificially and naturally contaminated food samples (Auvray *et al.* personal communication).

The last two systems may be also used as an extraction control by adding 100 copies of the pUC 19 plasmid or of the pIAC-STEC to the sample aliquot prior to the DNA purification step.

ANNEX E

The Real-time PCR protocol described is based on the use of the following primers and probes which shall be considered as reference reagents. However, other primers and probes may be used provided that they have been recognised equivalent to those indicated in the tables E.1 and E.2 according to the ISO 16140 rules.

Primers and probes for the PCR assays

Tables E.1 and E.2 provides respectively the primers and probes sequences for:

- the detection of *stx* and *eae* genes by real-time-PCR (PCR A);
- the detection of serogroup-related genes using real-time-PCR (PCR B).

In these tables, the chemistry of the reporter and quencher fluorophores is not indicated, being largely dependent on the real time PCR instruments available in each laboratory.

Table E.1: Degenerate primers and TaqMan probes used for 5' nuclease PCR assays. (§3 and *2)

Target gene	Forward primer, reverse primer and probe sequences (5'-3') ^a	Amplicon size (bp)	Location within sequence	GenBank accession number
<i>stx1</i> [§]	TTTGTACTGTSACAGCWGAAGCYTTACG CCCCAGTTCARWGTRAGRTCMACRTC Probe -CTGGATGATCTCAGTGGGCGTTCTTATGTAA	131	878–906 983–1008 941–971	M16625
<i>stx2</i> ^{§b}	TTTGTACTGTSACAGCWGAAGCYTTACG CCCCAGTTCARWGTRAGRTCMACRTC Probe -TCGTCAGGCACTGTCTGAAACTGCTCC	128	785–813 785–813 838–864	X07865
<i>eae</i> [*]	CAT TGA TCA GGA TTT TTC TGG TGA TA CTC ATG CGG AAA TAG CCG TTA Probe -ATAGTCTCGCCAGTATTCGCCACCAATACC	102	899-924 1000-979 966-936	Z11541

^a In the sequence Y is (C, T), S is (C, G), W is (A, T), R is (A, G), M is (A, C).

^b This combination of primer/probe recognises all the *stx2* variants but the *stx2f*

Table E.2. Primers and probes used for amplification of O antigen specific genes in 5' nuclease PCR assays. (§3 and * 4)

Target gene (serogroup)	Forward primer, reverse primer and probe sequences (5'-3')	Amplicon size (bp)	Location within sequence	GenBank accession number
[§] <i>rfbE</i> (O157)	TTTCACACTTATTGGATGGTCTCAA CGATGAGTTTATCTGCAAGGTGAT Probe -AGGACCGCAGAGGAAAGAGAGGAATTAAGG	88	348–372 412–435 381–410	AF163329
[§] <i>wbdI</i> (O111)	CGAGGCAACACATTATATAGTGCTTT TTTTTGAATAGTTATGAACATCTTGTTTAGC Probe -TTGAATCTCCCAGATGATCAACATCGTGAA	146	3464–3489 3579–3609 3519–3548	AF078736
[§] <i>wzx</i> (O26)	CGCGACGGCAGAGAAAATT AGCAGGCTTTTATATTCTCCAACCTT Probe -CCCCGTTAAATCAATACTATTTACGAGGTTGA	135	5648–5666 5757–5782 5692–5724	AF529080
[§] <i>ihp1</i> (O145)	CGATAATATTTACCCACCCAGTACAG GCCGCCGCAATGCTT Probe -CCGCCATTGAGAAATGCACACAATATCG	132	1383–1408 1500–1514 1472–1498	AF531429
[*] <i>wzx</i> (O103)	CAAGGTGATTACGAAAATGCATGT GAAAAAAGCACCCCCGTAATTAT Probe -CATAGCCTGTTGTTTTAT	99	4299–4323 4397–4375 4356–4373	AY532664

Annex F

Isolation of STEC strains

Follow the procedure described below to isolate STEC strains from real time PCR positive samples:

- 1) Perform a serogroup-specific enrichment (SSE) on the remaining enrichment culture (see Note 1)
- 2) Streak SSE onto TBX or other suitable medium (see note 2). Incubate for 18 to 24 hours at 37°C
- 3) Pick up 10 to 50 colonies with *E. coli* morphology or with characteristic aspect (see Note 5) and point-inoculate on nutrient agar (NA) (see Note 3) and H₂O (the colonies may be pooled in water up to a number of ten per pool).
- 4) Perform the detection of the *stx*-coding gene and the *eae* gene on the isolated colonies or the H₂O pools (see Note 4).
- 5) If a pool is positive, go back to NA and assay the individual colonies forming the positive pool in order to select one single positive colony.
- 6) Identify the colonies as *E. coli* and confirm the serogroup the sample was positive to in the screening PCR assay (e.g. by PCR B in the Annex E), see Note 5.
- 7) Isolates may be sent to the a Reference Laboratory for further characterization.

NOTE 1: Serogroup-specific enrichment may be achieved by using immunocapture systems such as immuno-magnetic separation (IMS) or equivalent. Generally, refer to the instruction supplied by the manufacturer.

For O157 positive samples, use ISO 16654 or alternative methods validated according to ISO 16140.

NOTE 2: For O157 positive samples, use ISO 16654 or alternative methods validated according to ISO 16140. Sorbitol-fermenting *E. coli* O157 are susceptible to tellurite contained in the CT SMAC medium indicated in ISO 16654. Therefore the use of a second SMAC isolation plate without antibiotics is recommended. In the absence of Sorbitol-negative colonies on the plates, the screening of Sorbitol-positive colonies is suggested.

For STEC O26 isolation, a differential solid media (MacConkey) containing Rhamnose instead of lactose is commercially available (RMAC). It is very effective in distinguishing STEC O26 strains, which do not ferment Rhamnose, from other *E. coli*.

NOTE 3: There are several types of nutrient agar media commercially available either ready to use plates or prepared in house from dehydrated powders. Every type of non-

selective nutrient agar media (e.g. TSA) is suitable for the purpose of maintaining the colonies for further characterisation. Enterohaemolysin Agar, can also be used. It gives the advantage to detect the Enterohaemolysin production, which is a common feature of STEC pathogenic to humans.

NOTE 4: The Real Time PCR described in this protocol may be adopted to confirm the presence of the *stx* and *eae* in the isolated strains. Conventional PCR may be used as an alternative (<http://www.iss.it/vtec/work/cont.php?id=152&lang=2&tipo=3>).

NOTE 5: Colony confirmation as *E. coli* may be achieved by using any commercial biochemical multi-assay or by assessing the indole production. Confirmation of the serogroup may be achieved either by PCR or by agglutination with commercial antisera.

Bibliography

- [1] Monitoring of verotoxigenic *Escherichia coli* (VTEC) and identification of human pathogenic VTEC types Scientific Opinion of the Panel on Biological Hazards (Question No EFSA-Q-2007-036) *The EFSA Journal* (2007) 579, 1-61
- [2] Nielsen EM, Andersen MT. Detection and characterization of verocytotoxin-producing *Escherichia coli* by automated 5' nuclease PCR assay. 2003 *J Clin Microbiol.* 41(7):2884-2893.
- [3] Perelle S, Dilasser F, Grout J, Fach P. Detection by 5'-nuclease PCR of Shiga-toxin producing *Escherichia coli* O26, O55, O91, O103, O111, O113, O145 and O157:H7, associated with the world's most frequent clinical cases. 2004 *Mol Cell Probes.* 18(3):185-192.
- [4] Perelle S, Dilasser F, Grout J, Fach P. Detection of *Escherichia coli* serogroup O103 by real-time polymerase chain reaction. 2005 *J Appl Microbiol* 98(5):1162-1168.
- [5] Vimont A, Delignette-Muller ML, Vernozy-Rozand C. Supplementation of enrichment broths by novobiocin for detecting Shiga toxin-producing *Escherichia coli* from food: a controversial use. *Lett Appl Microbiol* 2007 44(3):326-331.
- [6] Uemura R, Sueyoshi M, Nagayoshi M, Nagatomo H. Antimicrobial susceptibilities of Shiga toxin-producing *Escherichia coli* isolates from pigs with edema disease in Japan. 2003 *Microbiol Immunol* 47(1):57-61.
- [7] Sambrook and Russel, *Molecular cloning, a laboratory manual*. Cold Spring Harbor Laboratory Press. third edition, 2001

- [8] Paton AW, Paton JC. Detection and characterization of Shiga *toxigenic Escherichia coli* by using multiplex PCR assays for *stx1*, *stx2*, *eaeA*, enterohemorrhagic *E. coli hlyA*, *rfbO111*, and *rfbO157*. 1998 *J Clin Microbiol* 36: 598-602.
- [9] O'Brien AD, Newland JW, Miller SF, Holmes RK, Smith HW, Formal SB. Shiga-like toxin-converting phages from *Escherichia coli* strains that cause hemorrhagic colitis or infantile diarrhea. 1984 *Science* 226: 694-696.
- [10] Perna N.T. *et al.* Genome sequence of enterohaemorrhagic *Escherichia coli* O157:H7. 2001 *Nature* 409: 529-533.
- [11] Fricker M, Messelhäusser U, Busch U, Scherer S, Ehling-Schulz M. Diagnostic real-time PCR assays for the detection of emetic *Bacillus cereus* strains in foods and recent food-borne outbreaks. 2007 *Appl Environ Microbiol* 73(6):1892-1898.
- [12] Hoorfar J, Malorny B, Abdulmawjood A, Cook N, Wagner M, Fach P. Practical considerations in design of internal amplification controls for diagnostic PCR assays. 2004 *J Clin Microbiol* 42(5):1863-1868.
- [13] Beutin L, Jahn S, Fach P. Evaluation of the 'GeneDisc' real-time PCR system for detection of enterohaemorrhagic *Escherichia coli* (EHEC) O26, O103, O111, O145 and O157 strains according to their virulence markers and their O- and H-antigen-associated genes. 2009. *J Appl Microbiol.* 106:1122-32.